


VITROCELL® Cloud α PowerVent

Increase of nebulizations and evacuation of residual gaseous compounds


A typical Cloud exposure allows 3 consecutive nebulizations. The PowerVent option enables an increase of the dose by adding further nebulizations. After each nebulization the humidity is removed by a short venting period before starting the new one. The venting function can be edited in the software. The PowerVent unit is delivered ready-to-use with an integrated vacuum pump.

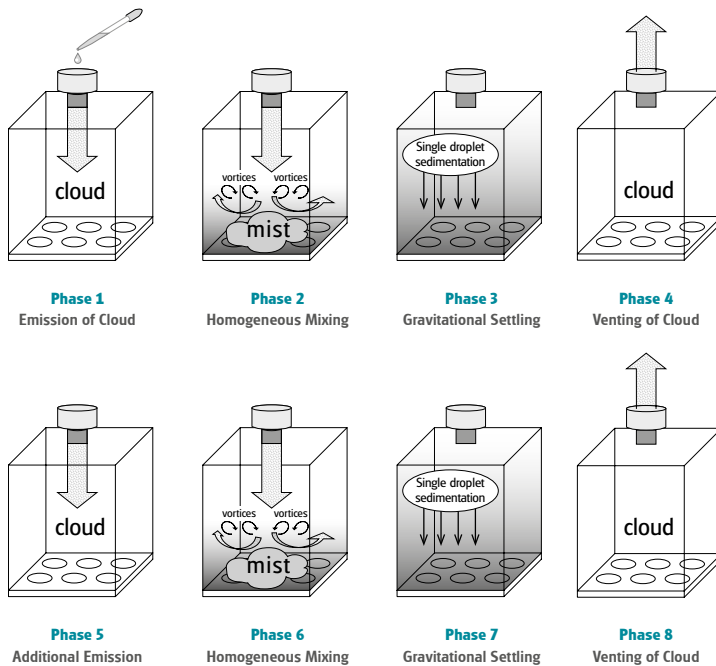

Vitrocell® Cloud Alpha 6 with PowerVent option.

PowerVent components

- (1) HEPA filter for filtered air
- (2) PowerVent adapter
- (3) HEPA filter
- (4) PowerVent Unit with vacuum pump

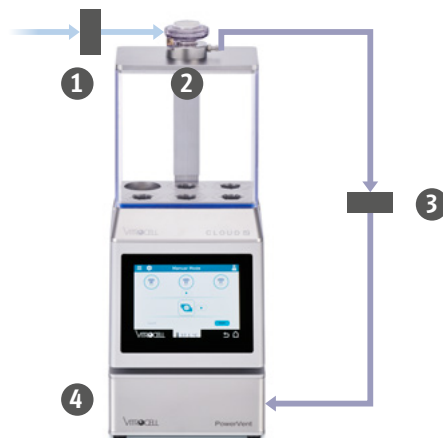

Key Features:

- Available for all Cloud Alpha systems
- Fast removal of residual gaseous compounds
- Allows >3 consecutive nebulizations
- Integrated vacuum pump
- Function can be edited in the software